



Resorbable Haemostatic Sponge

Absorbable Gelatin Sponge USP

HEMOSPONGE®



ISO 13485:2016
CERTIFIED COMPANY



INDIAN FDA
APPROVED



e- Catalogue

Power of 3D MATRIX



FAST AND EFFECTIVE HAEMOSTAT

Stops oozing blood in 2-4 minutes




ACTIVATES NATURAL CLOTTING PROCESS ONLY

Provides surface to auto activate platelets



HIGHLY ABSORBENT

Absorbs and holds the blood upto 45 times the weight of sponge

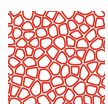

HIGH SURFACE AREA AND UNIFORM POROSITY MATRIX

Faster clotting and better clot holding capacity



100% RESORBABLE

Completely resorbed in body within four weeks



PLIABLE

Product can be rolled, cut, shaped to fit bleeding site



BIO-COMPATIBLE

Kinder to body tissues (complies ISO 10993 – Biocompatibility Standard)



PURE AND SAFE

Made up of pure pharma grade gelatin only, no chemicals, no drugs



pH NEUTRAL

Can be soaked with thrombin or other drugs for advanced effect

**USP 42
NF 37**

ASSURING QUALITY


Product complies US Pharmacopoeia 42, NF 37 quality standard

Product VARIANTS

	Variant	Sponge Size (mm)	Product reference	Selling Pack sizes
Regular Sizes	STANDARD	80 x 50 x 10	GL851	2 Pcs/Pack, 10 Pcs/Pack 12 Pcs/Pack, 20 Pcs/Pack
	EXPORT	70 x 50 x 10	GL751	2 Pcs/Pack, 10 Pcs/Pack 12 Pcs/Pack, 20 Pcs/Pack
Dental Range	DENTAL SL	10 x 10 x 10	GL111	10 Pcs/Pack, 24 Pcs/Pack, 32 Pcs/Pack, 50 Pcs/Pack
	DENTAL M	14 x 14 x 07	GL1414	10 Pcs/Pack, 50 Pcs/Pack
	DENTAL L	15 x 15 x 15	GL1515	10 Pcs/Pack, 50 Pcs/Pack
	DENTAL XL	20 x 20 x 07	GL202007-D	20 Pcs/Pack, 48 Pcs/Pack
	DENTAL MAX	30 x 30 x 10	GL303010-D	10 Pcs/Pack, 20 Pcs/Pack
ENT Range	ENT SL	20 x 20 x 07	GL202007-E	20 Pcs/Pack, 48 Pcs/Pack
	ENT XL	30 x 30 x 10	GL303010-E	10 Pcs/Pack, 20 Pcs/Pack
	EPISTAXIS	80 x 15 x 08	GL10015	10 Pcs/Pack
	NASAL TAMPON	80 x 10 (Dia)	GL8010	10 Pcs/Pack, 20 Pcs/Pack
	WICK SL	15 x 09	GL1509	10 Pcs/Pack, 20 Pcs/Pack
	WICK XL	24 x 09	GL2409	10 Pcs/Pack, 20 Pcs/Pack
Tampon	TAMPON	80 x 30 (Dia.)	GL8030	2 Pcs/Pack, 6 Pcs/Pack, 10 Pcs/Pack
	TAMPON – SL	60 x 30 (Dia.)	GL6030	2 Pcs/Pack, 6 Pcs/Pack, 10 Pcs/Pack
Powder	POWDER	1 Gram	GL001	6 Pcs/Pack
Film Range	SPECIAL – 1	80 x 50 x 01	GL805001	10 Pcs/Pack
	SPECIAL – 2	70 x 50 x 01	GL705001	10 Pcs/Pack
	FILM – 80	80 x 25 x 01	GL80	12 Pcs/Pack
	FILM – 120	120 x 100 x 01	GL120	12 Pcs/Pack
	FILM – 240	240 x 120 x 01	GL240	12 Pcs/Pack
	FILM – 200	200 x 70 x 0.5	GL200	12 Pcs/Pack

Custom Size : HEMOSPONGE is also available in many sizes other than mentioned above, please contact on info@goodwill-lifesciences.com for more details.

Hemosponge REGULAR SIZES

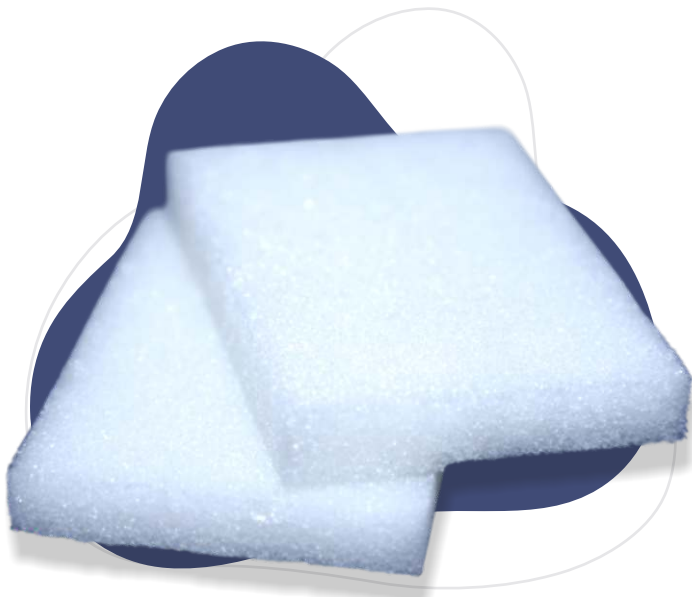

- Rectangle shaped sponge with uniform porosity
- Available in wide variety of packing size
- Can be cut, rolled or mould to desired shape to fit any bleeding site

**General surgeries | Laparoscopic Surgeries | Neurology
Cardiology | Orthopaedic | Nephrology | Oncological
Surgeries | Gynaecology | Proctology | Oral, Dental and
ENT**



**For all Surgical
procedures**



**Cut, Rolled,
Shaped**

Hemosponge TAMPONS



- Cylindrically-shaped sponge with hollow region in centre
- Easily inserted with the help of proctoscope
- Absorbs continuously oozing fluids for longer period
- Enhances patient comfort by reducing post-operative pain and irritation
- Expelled easily with first natural bowl movement



**Coloproctological
Surgeries**




**Post -
Hemorrhoid
pain relief**



**Gynecological
Surgeries**

Hemosponge DENTAL RANGE

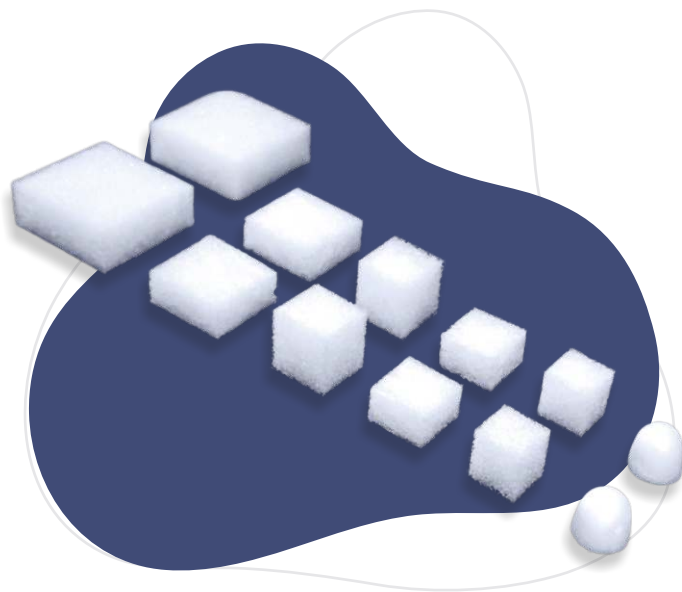

- Small size sponges for use in simple extraction to orthodontic and maxillofacial -oral surgery
- No need to cut – easily fit bleeding sockets
- Better than cotton based products – No chance of re-bleeding
- Can be soaked with thrombin for quick haemostasis
- Can be soaked with antibiotic – for local antiseptic effect
- No need to remove after surgery – provides higher patient comfort
- Individual and sterile piece packing – ready to use and cost – effective
- Made up of pure gelatin – Harmless if swallowed



**Cost
Effective**



**For All Dental
Procedures**

**SMALL
BIG**

**Small Size
Big Effect**



**Replacement
of Cotton**

Hemosponge ENT RANGE




- Specially shaped sponges to fit into small surgical cavities
- Widely accepted and effective haemostatic – packing material
- 100% absorbable – avoids patient discomfort associated with removal of non-absorbable packing material.
- Works as a scaffold to support growth of near-by tissue and fastens the healing process of ENT surgeries.



**Special Sizes to
fit ENT cavities**

Hemosponge **FILM RANGE**

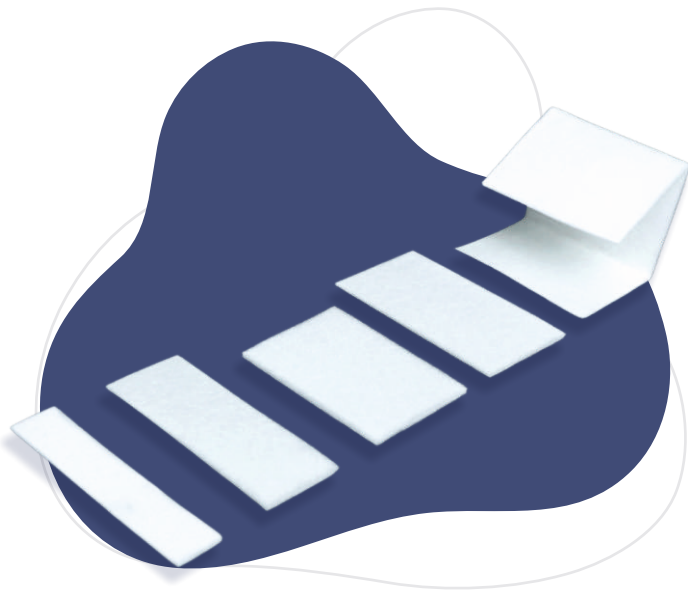

- Useful for surgeries where large size of sponge is impossible to place i.e. in laminectomy
- Designed for use in thoracic, ocular and neurosurgery
- Available in various sizes with thickness ranging from 0.5 to 1.00 mm

Hemosponge **POWDER**



- Highly effective hemostatic for the sites with uneven surface
- Mixed with saline and kneaded to form putty of desired consistency before application
- The putty has hemostatic effect and is applied directly to bleeding surface



**For Uneven
Surface**

About COMPANY

- An ISO 14644 Cleanrooms Production Line
- An ISO 13485:2016 Certified QMS
- CDSCO Approved Class D Manufacturing Unit
- In House Quality Control Laboratories



E- Catalogue



All variants



Dental Range



ENT Range



Follow Us For
Updates

Get in **TOUCH**



www.goodwill-lifesciences.com



info@goodwill-lifesciences.com | sales@goodwill-lifesciences.com



+91-87586 18110 | +91-97147 44173



Goodwill Lifesciences Private Limited

283/284/285, Block no. 392/A, Anjani Textile Park - 4, Segwa-Paria Road, At- Paria - 394130,
Ta. Olpad, Surat, Gujarat, India, Mfg Lic no.: MFG/MD/2022/000029